


LACROSSE
ApertaTM_{NSE}
Coronary Dilatation Catheter



Catéter de
Dilatación Coronaria

 **NIPRO**



La nueva generación de NSE

Un diseño modificado de NSE scoring balloon garantiza un mayor efecto scoring, en una plataforma no compatible, con un perfil de cruce más bajo.

CARACTERÍSTICAS

Aperta NSE

Elementos : integrados

NSE ALPHA

Elementos : libres - flotando



Elementos Integrados

La integración de los elementos con el balón mantiene la posición del balón con una separación de 120 grados durante la expansión del globo.

Alta Presión

Una plataforma que no cumple proporciona un aumento en la transmisión de empuje desde el vértice del elemento.


NP

14atm

RBP

24atm

NSE ALPHA

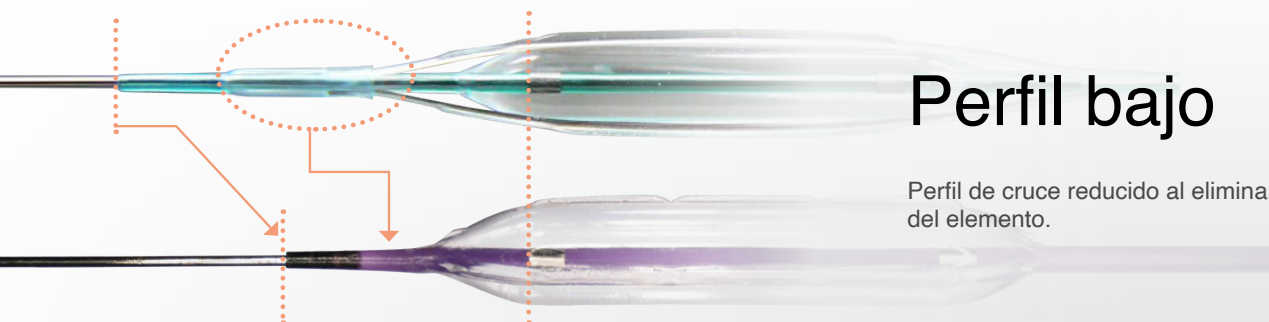

Aperta NSE

Perfil bajo

Perfil de cruce reducido al eliminar la conexión del elemento.

DISEÑO DEL CATÉTER

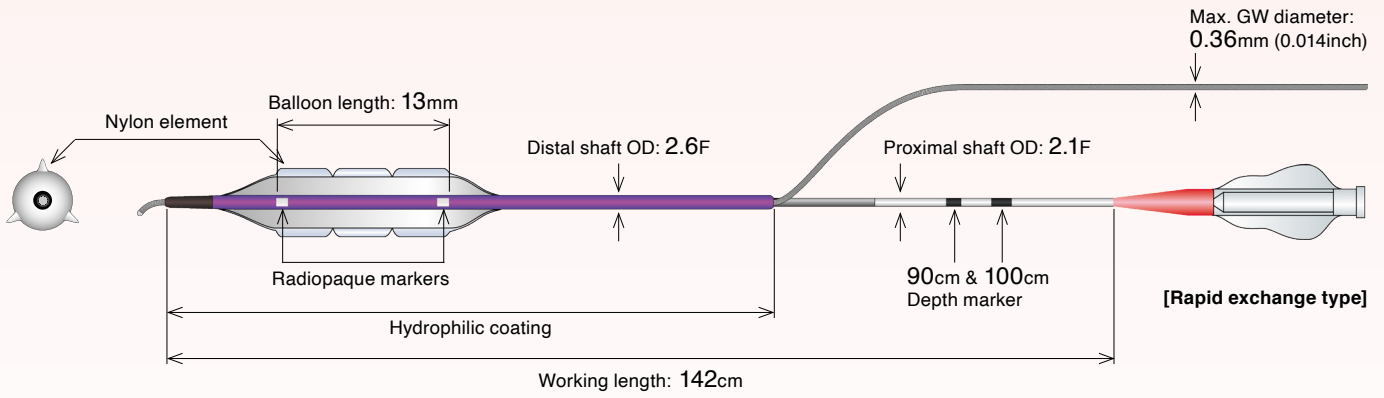


TABLA DE CUMPLIMIENTOS

Inflation pressure		Balloon diameter (mm)							
atm	kPa	2.00mm	2.25mm	2.50mm	2.75mm	3.00mm	3.25mm	3.50mm	4.00mm
4	4x10 ²	1.78	2.02	2.27	2.45	2.70	2.95	3.18	3.50
5	5x10 ²	1.81	2.05	2.30	2.50	2.73	2.98	3.23	3.55
6	6x10 ²	1.84	2.08	2.33	2.55	2.76	3.01	3.26	3.60
7	7x10 ²	1.86	2.11	2.36	2.59	2.79	3.04	3.29	3.65
8	8x10 ²	1.88	2.13	2.38	2.63	2.82	3.07	3.32	3.70
9	9x10 ²	1.90	2.15	2.40	2.65	2.85	3.10	3.35	3.75
10	10x10 ²	1.92	2.17	2.42	2.67	2.88	3.13	3.38	3.80
11	11x10 ²	1.94	2.19	2.44	2.69	2.91	3.16	3.41	3.85
12	12x10 ²	1.96	2.21	2.46	2.71	2.94	3.19	3.44	3.90
13	13x10 ²	1.98	2.23	2.48	2.73	2.97	3.22	3.47	3.95
14 (NP)	14x10²	2.00	2.25	2.50	2.75	3.00	3.25	3.50	4.00
15	15x10 ²	2.02	2.27	2.52	2.77	3.02	3.27	3.52	4.02
16	16x10 ²	2.04	2.29	2.54	2.79	3.04	3.29	3.54	4.04
17	17x10 ²	2.06	2.31	2.56	2.81	3.06	3.31	3.56	4.06
18	18x10 ²	2.08	2.33	2.58	2.83	3.08	3.33	3.58	4.08
19	19x10 ²	2.09	2.35	2.60	2.85	3.10	3.35	3.60	4.10
20	20x10 ²	2.10	2.36	2.62	2.87	3.12	3.37	3.62	4.12
21	21x10 ²	2.11	2.37	2.64	2.89	3.14	3.39	3.64	4.14
22	22x10 ²	2.12	2.38	2.66	2.91	3.16	3.41	3.66	4.16
23	23x10 ²	2.13	2.39	2.68	2.93	3.18	3.43	3.68	4.18
24 (RBP)	24x10²	2.14	2.40	2.70	2.95	3.20	3.45	3.70	4.20
25	25x10 ²	2.15	2.41	2.72	2.97	3.22	3.47	3.72	4.22
26	26x10 ²	2.16	2.42	2.74	2.99	3.24	3.49	3.74	4.24
27	27x10 ²	2.17	2.43	2.76	3.01	3.26	3.51	3.76	4.26

INFORMACIÓN PARA ORDENAR

Model no.	Balloon diameter (mm)	Balloon length (mm)	Proximal shaft OD (F / mm)	Distal shaft OD (F / mm)	Nylon elements
APN20013	2.00	13	2.1 / 0.68	2.6 / 0.86	3
APN22513	2.25				
APN25013	2.50				
APN27513	2.75				
APN30013	3.00				
APN32513	3.25				
APN35013	3.50				
APN40013	4.00				



https://med.nipro.co.jp/en_products

Aperta NSE Coronary Dilatation Catheter
Manufacturer: Goodman Co., Ltd.

- ▶ Please consult your sales representative regarding stock availability.
- ▶ Lacrosse and Aperta NSE are trademarks of Goodman Co., Ltd.
- ▶ Please refer to the Instructions For Use for further details regarding product use and specifications.